
Slides for School Model



Test No. 1, SF-1 — Binocularity Both Eyes

To determine whether both eyes are being used together.

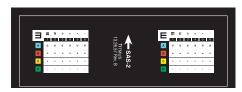

Test No. 2, SAS-2 — Acuity Right and Left Eye

Conventional Tumbling "E"symbol.
Child states which way the "E"points in each square.
Acuity levels of 20/20, 20/30, 20/40 and 20/50.



Test No. 3, SPS-3 — Heterophorias

Muscle imbalance – vertical and lateral – is tested at predetermined pass/fail levels.

This test is administered first at far, then at near.



Test No. 4, SPCD-1 - Color Perception

Titmus Pseudo-Isochromatic plates reliably detect gross color anomalies, but do not classify the type of deficiency.



Test No. 5, SA-BRL-4 — Visual Acuity

Results calibrated at 20/20, 20/30, 20/40, 20/50, 20/70, 20/100, and 20/200.

Tests can be administered binocularly or with either eye occluded.



Test No. 6, SLP-1 — Lateral Phoria

Measures Esophoria or Exophoria in 1 prism diopter steps.



Test No. 7, SVP-1 — Vertical Phoria

Measures right or left Hyperphorias in one-half prism diopter steps.



Test No. 8, SA-BRL-2A — Visual Acuity

Results calibrated at 20/20, 20/30, 20/40, 20/50, 20/70 and 20/100.

Tests can be administered binocularly or

Tests can be administered binocularly or with either eye occluded.

Improvements have been made to test slides to enhance viewing performance. Patent Pending.



