
Slides for Professional Model



Test No. 1, SA-APS-1 Allen Picture Test

Line 1 – Binocularity – Tests performance of both eyes together. Lines 2, 3, 4 – Acuity – Right and Left Eye Calibrated at 20/30, 20/40 and 20/50.

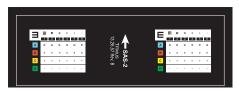

Test No. 2, SAS-2 — Acuity Right and Left Eye

Conventional Tumbling "E"symbol.
Child states which way the "E"points in each square.
Acuity levels of 20/20, 20/30, 20/40 and 20/50.



Test No. 3, SPS-3 — Heterophorias

Muscle imbalance – vertical and lateral – is tested at predetermined pass/fail levels.

This test is administered first at far, then at near.



Test No. 4, SCD-1 — Color Perception

Consists of six accurately and authentically reproduced Ishihara Pseudo-Isochromatic plates. Detects color deficiency but does not classify.



Test No. 5, SA-BRL-4 — Visual Acuity

Results calibrated at 20/20, 20/30, 20/40, 20/50, 20/70, 20/100, and 20/200.

Tests can be administered binocularly or with either eye occluded.



Test No. 6, SLP-1 - Lateral Phoria

Measures Esophoria or Exophoria in 1 prism diopter steps.



Test No. 7, SVP-1 - Vertical Phoria

Measures right or left Hyperphorias in one-half prism diopter steps.



Test No. 8, SA-BRL-2A — Visual Acuity

Results calibrated at 20/20, 20/30, 20/40, 20/50, 20/70, and 20/100.

Tests can be administered binocularly or with either eye occluded.

Improvements have been made to test slides to enhance viewing performance. Patent Pending.



